
GESKE

GERMAN BEAUTY TECH

Quickstart Guide

www.geske.com

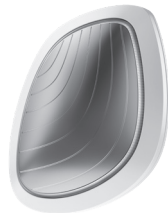
GESKE Beauty Tech GmbH,
Leipziger Platz 18, 10117 Berlin, Germany
info@geske.com, www.geske.com

SmartAppGuided™
Sonic Warm & Cool Mask | 9 in 1

Short Introduction


Scan for the full instructions
in your language

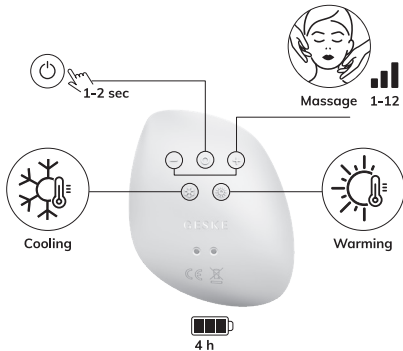GESKE German Beauty Tech App
Available in 36 Languages

Apple and the Apple logo are trademarks of Apple Inc.
Google Play and the Google Play logo are trademarks of Google LLC.

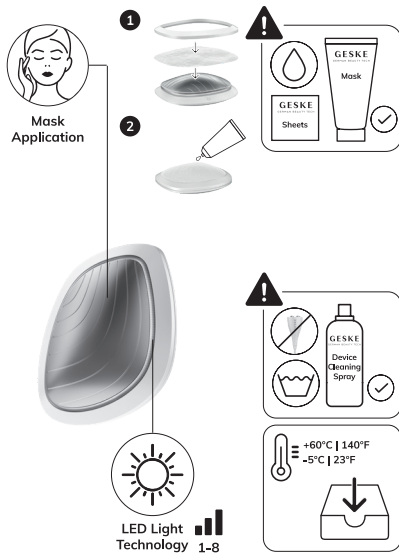**Read the safety information before first use.
Keep both instructions for use.**

YOUR PRODUCT AT A GLANCE

Your product is a multi-functional mask applicator. It is safe to be used up to 3 times a day, please do not use more frequently. Please do not use the light and pulsation mode more than once per day. For maximal safety and ideal results, always follow your personalized Skin Care Routine created by the GESKE App.

Download the GESKE App and get your personalized routine.

SPECIFICATIONS:

Materials: ABS, silicone, PC, aluminium alloy | **Storage temperature:** -5 - +60 °C/23 - 140 °F | **Battery:** Li-Ion 3,7 V | **Battery Capacity:** 800 mAh | **Maximum Charging Current:** 450 mA | **Charging Voltage:** 5V | **Charging Time:** approx. 4 hours | **Operating Time:** Operating time differs depending on activated technologies & intensities. | **Interface:** USB | **Programs:** Vibration: 12 intensities, Heating, Cooling, LED Light Technology: 8 colors

IMPORTANT SAFETY NOTES! READ BEFORE USAGE!

The product is not suitable for use under the following conditions:

- Do not use if you are known to be hypersensitive to heat (e.g. heat urticaria)
- If a feeling of (too) intense heat occurs during usage, it must be stopped immediately
- The safety and effectiveness of the product has never been assessed in conjunction with non-GESKE skin care products and masks. We recommend using the product only with GESKE products. Use with other products at your own risk
- The product can emit intense light, avoid applying it in the direction of or near an open eye
- Do not use the product on irritated or injured skin or open wounds (after any kind of injury), in the case of bacterial or viral infections of the skin (either in the usage area and elsewhere)
- This product is not intended to treat medical skin conditions or diseases. If you have medical skin conditions or diseases, consult your dermatologist before using this product
- Do not use the product if it is 1) untypically warm, 2) mechanically damaged or deformed, 3) discolored
- If pain and/or discomfort occurs during use, discontinue use and seek medical advice
- The eye area is very sensitive: Be extra careful in this area and avoid any contact of the product with the eyelids or the eye itself

- Only use the product as recommended in your personalized routine
- For hygiene reasons, we recommend you do not share your product with anyone else
- The product is intended for cosmetic use on the face. GESKE assumes no liability for harmful consequences arising from the use on other body areas or other improper use, connection to unsuitable voltage sources, contaminated conducting liquid or balls
- Always make sure the plug and socket are completely dry before charging. Otherwise, electric shocks, a short circuit or fire may result. The product must not be used while connected to the power outlet

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