GESKE

GERMAN BEAUTY TECH

Quickstart Guide

SmartAppGuided™
Sonic Thermo Facial
Brush & Face-Lifter | 8 in 1

www.geske.com

GESKE Beauty Tech GmbH,
Leipziger Platz 18, 10117 Berlin, Germany
info@geske.com, www.geske.com

Short Introduction

Scan for the full instructions
in your language

Massage 1-18

1-2 sec

Warming


Warming

MicroCurrent 1-5

+60°C | 140°F
-5°C | 23°F

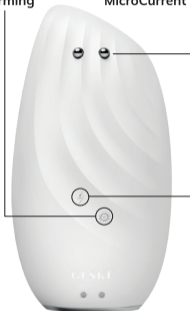4.0 h

GESKE German Beauty Tech App
Available in 36 Languages

Apple, the Apple logo and Apple Watch are trademarks of Apple Inc.
Google Play and the Google Play logo are trademarks of Google LLC.

The safety information must also be read before first use. Keep both instructions for use.

YOUR PRODUCT AT A GLANCE

Your product is a multi-functional Facial Brush. It is safe to be used up to 3 times a day, please do not use more frequently. Please do not use the MicroCurrent-Lifting modes more than once per day. For maximal safety and ideal results, always follow your personalized Skin Care Routine created by the GESKE App.

Download the GESKE App and get your personalized routine.

Intended Use: The Geske Sonic Thermo Facial Brush & Face-Lifter Device is an over-the-counter home use device, which is intended for thermo facial cleansing of dirt, oil, & dead skin, and facial skin stimulation with MicroCurrent technology

SPECIFICATIONS:

Material: ABS, silicone, aluminium alloy, brass | **Storage temperature:** -5 - +60°C | **Battery:** Li-Ion 3.7V | **Battery Capacity:** 1200 mAh | **Maximum Charging Current:** 440 mA | **Charging Voltage:** 5V | **Charging Time:** approx. 4,0 hours | **Operating Time:** Operating time differs depending on activated technologies & intensities | **Interface:** USB | **Programs:** Vibration: 18 intensities, Heating, MicroCurrent: 5 intensities

IMPORTANT SAFETY INFORMATION! READ BEFORE USE!

The product is not suitable for use under the following conditions:

- Do not use if you are known to be hypersensitive to heat (e.g. heat urticaria).
- If a feeling of (too) intense heat occurs during the treatment, the treatment must be stopped immediately.
- If you have an implanted medical device or other electronic device such as a pacemaker
- After laser treatment or chemical peeling
- After cosmetic facial surgery

In the following cases, a consultation with a physician is necessary as to whether the product is suitable for you to use:

- If you have a heart disease
- In case of pregnancy. The safety of use during pregnancy has not been established

Further information on use:

- Before using the Microcurrent Facelift Technology, the product must be completely dry
- Use of the Microcurrent Facelift Technology only in combination with the GESKE serum

- Do not use the Microcurrent function on the chest area, around the eyes, midline of the neck, genital/groin area, or on red and inflamed skin, over birthmarks and enlarged blood vessels
- Nothing is known about the long-term effects of Microcurrent
- As a result of stimulation of the optic nerve, a flashing light may be perceived during use. If this persists when not using the product, seek medical advice
- A slight tingling sensation may occur during application. Decrease the intensity of the Microcurrent to reduce or eliminate tingling
- Electronic monitoring devices may not function properly while the device is in use
- Never use the product while driving or operating heavy machinery
- Do not poke any sharp objects into the product's openings
- The product may only be used with a safety extra-low voltage (SELV) adapter
- We recommend using standard IEC60950 power supplies to charge the device
- Do not use the product on irritated or injured skin or open wounds (after any kind of injury), in the case of bacterial or viral infections of the skin (either in the treatment area and elsewhere)
- Only use the product if you have skin problems and diseases or other medical concerns after consulting a dermatologist
- Do not use the product if it is 1) damaged or deformed, 2) discolored.
- If pain and/or discomfort occurs during use, discontinue use and seek medical advice
- The eye area is very sensitive: Be extra careful in this area and avoid any contact of the product with the eyelids or the eye itself
- Only use the product as recommended in your personalized routine
- For hygiene reasons, we recommend you do not share your product with anyone else
- The product is intended for cosmetic use on the face. GESKE assumes no liability for harmful consequences arising from the use on other body areas or other improper use, connection to unsuitable voltage sources, contaminated conducting liquid or balls
- Always make sure the plug and socket are completely dry before charging. Otherwise, electric shocks, a short circuit or fire may result. The product must not be used while connected to the power outlet

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