# GESKE

GERMAN BEAUTY TECH

## Quickstart Guide

[www.geske.com](http://www.geske.com)

GESKE Beauty Tech GmbH,  
Leipziger Platz 18, 10117 Berlin, Germany  
[info@geske.com](mailto:info@geske.com), [www.geske.com](http://www.geske.com)

SmartAppGuided™  
**MicroNeedle Face & Body Roller | 9 in 1**

# Short Introduction

Scan for the full instructions  
in your language



Massage

2 sec

before  
use

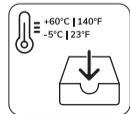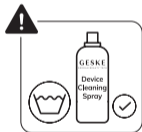GESKE German Beauty Tech App  
Available in 45 Languages

Apple and the Apple logo are trademarks of Apple Inc.  
Google Play and the Google Play logo are trademarks of Google LLC.

EN

**Read the safety information before first use.**

**Keep both instructions for use.**

## YOUR PRODUCT AT A GLANCE

Your product is a multi-functional MicroNeedle roller. It is safe to be used up to 2 times per week, at maximum 3 times per week. Please do not use more frequently. This limitation only applies to the most intense function of the device: microneedling. For maximal safety and ideal results, always follow your personalized Skin Care Routine created by the GESKE App.

**Download the GESKE app and get your personalized routine.**

### SPECIFICATIONS:

**Material:** silicone, ABS, PC, stainless steel, rose quartz |  
**Storage temperature:** -5 - +60°C/23 - 140°F | **Battery:** Li-Ion 3,7V |  
**Battery Capacity:** 140 mAh | **Maximum Charging Current:** 140 mA |  
**Charging Voltage:** 5V | **Charging Time:** approx. 2,5 hours |  
**Operating Time:** Operating time differs depending on activated technologies & intensities. | **Interface:** USB | **Programs:** Vibration: 1 speed level, LED Light Technology: 1 colour

### IMPORTANT SAFETY NOTES! READ BEFORE USAGE!

**The product is not suitable for use under the following conditions:**

- If you suffer from skin cancer or precancerous stages (actinic keratoses)
  - If you suffer from hemorrhagic diseases (congenital blood clotting disorders) and/or take medication that affects blood clotting
  - The safety and effectiveness of the product has never been assessed in conjunction with non-GESKE skin care products and masks. We recommend using the product only with GESKE products. Use with other products at your own risk
  - After laser treatment or chemical peeling
  - After facial surgery, both cosmetically and medically indicated
- In the following cases, a consultation with a physician is necessary as to whether the product is suitable for you to use:**
- If you have a heart disease
  - In case of pregnancy. The safety of use during pregnancy has not been established
  - If you have derma fillers or injections

### Further information on use:

- Only apply light pressure when using the microneedling head
- Microneedling can lead to mild swelling or bruising. Moisturize thoroughly after each session
- Avoid exercising, swimming and sun exposure for 24 h after a microneedling session. If you cannot avoid sun exposure apply a high SPF sunscreen.

Do not apply make-up for 24 h

- Avoid using products with acids and fragrance, as well as physical scrubs and exfoliators during and after a microneedling session
- Using products with active ingredients such as vitamin C or retinol may lead to more intense irritation of the skin. If this irritation persists for a longer period of time, seek medical advice
- Sanitize the needle attachment after each use with alcohol spray and let it air dry
- Do not use the product on irritated or injured skin or open wounds (after any kind of injury), in the case of bacterial or viral infections of the skin (either in the usage area and elsewhere)
- Do not use the product on nevi, moles or other (benign) skin lesions
- This product is not intended to treat medical skin conditions or diseases. If you have medical skin conditions or diseases, consult your dermatologist before using this product
- Do not use the product if it is 1) untypically warm, 2) mechanically damaged or deformed, 3) discolored
- If pain and/or discomfort occurs during use, discontinue use and seek medical advice
- The product can emit intense light, avoid applying it in the direction of or near an open eye
- The eye area is very sensitive: Be extra careful in this area and avoid any contact of the product with the eyelids or the eye itself
- Only use the product as recommended in your personalized routine
- For hygiene reasons, we recommend not sharing your product with others and replacing attachments as needed
- The product is intended for cosmetic use on the face. GESKE assumes no liability for harmful consequences arising from the use on other body areas or other improper use, connection to unsuitable voltage sources, contaminated conducting liquid or balls
- Always make sure the plug and socket are completely dry before charging. Otherwise, electric shocks, a short circuit or fire may result. The product must not be used while connected to the power outlet

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